

# Clear Control<sup>3D</sup>™

Sharpness is one thing, control is everything™

## Sharpness and Control of Incision

- 35% lower penetration force required than standard slit knives<sup>1</sup>
- Blade design provides tactile feel and incision control
- Ultra thin blade allows for less vertical tissue displacement for a more self-sealing wound

## 3D Incision


- Designed to facilitate instrument insertion through a multi-dimensional incision
- Micro bevels are designed to add vertical dimension, placing less stress on wound margins
- Proprietary Infinite Edge™ Technology produces unique designs and ultimate consistency

 INFINITE EDGE™  
TECHNOLOGY

 SharpPoint™  
The Surgeon's Edge™


# Penetration Force and Blade Thickness

Clear Control 3D™ Knife vs. Standard Slit Knife<sup>1</sup>

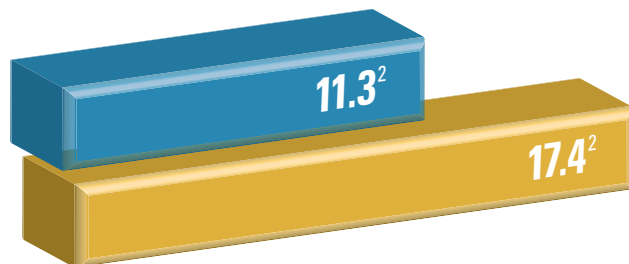

**35% Lower Penetration Force Required**

- Clear Control 3D™ blade
- Standard bevel-up slit blade

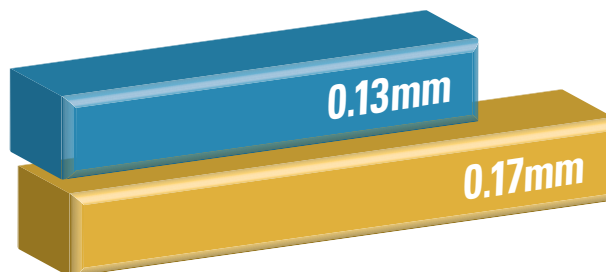


**22% Thinner**

<sup>1</sup>Comparative study on file of Surgical Specialties Corporation knives. Average values reported. <sup>2</sup>Single pass measured in grams.

## Clear Control 3D™ Slit Knives

SharpPoint™ Clear Control 3D™ slit knives are designed to facilitate instrument insertion through a multi-dimensional incision. Micro bevels are designed to add vertical dimension, placing less stress on wound margins, allowing for instrument insertion and manipulation. 6 per box.

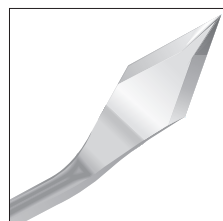

Knife	Blade/Tip	Description
<b>72-1550</b>	•	1.5mm Angled, full handle
<b>72-1650</b>	•	1.6mm Angled, full handle
<b>72-1850</b>	•	1.8mm Angled, full handle
<b>72-2050</b>	•	2.0mm Angled, full handle
<b>72-2250</b>	•	2.2mm Angled, full handle
<b>72-2450</b>	•	2.4mm Angled, full handle
<b>72-2550</b>	•	2.5mm Angled, full handle
<b>72-2650</b>	•	2.65mm Angled, full handle
<b>72-2750</b>	•	2.75mm Angled, full handle
<b>72-2850</b>	•	2.85mm Angled, full handle
<b>72-3050</b>	•	3.0mm Angled, full handle
<b>72-3250</b>	•	3.2mm Angled, full handle

To learn more about the Clear Control 3D™ Slit Knife:

visit

- [www.AngioEduPRO.com](http://www.AngioEduPRO.com), our educational website for healthcare professionals
- Or use your favorite QR code scanning program on your mobile device

call

- 1 877 991-1110 (USA and Canada only) or +001 732 626-6466 (International)

email

- [SharpPoint@angio.com](mailto:SharpPoint@angio.com)



Angiotech, 100 Dennis Drive, Reading, PA 19606, U.S.A.  
732 626-6466 • 877 991-1110 (USA and Canada only) • Fax: 610 404-4010 • [www.SharpPoint.com](http://www.SharpPoint.com) • OPG-PM-437R1 03/13  
SharpPoint™, The Surgeon's Edge™, Infinite Edge™, and Clear Control 3D™ are trademarks of Surgical Specialties Corporation.  
© 2013 Angiotech Pharmaceuticals, Inc. © 2013 Surgical Specialties Corporation. All Rights Reserved.

